



Translux® Wave

Power and speed at your fingertips.

Lightweight, cordless and convenient.

The Translux Wave is a new state-of-the-art curing light that delivers optimum polymerization for all common camphorquinone based dental materials. Its lightweight, cordless and ergonomic design provides a safe and easy method for consistent light-curing.


Giving a hand to oral health.



KULZER
MITSUI CHEMICALS GROUP

Translux® Wave

- **Pen-style body:** The lightweight and ergonomic pen-style design enables effortless handling even for small hands. The advanced angle of the 360° rotatable light guide allows better access to the molar regions.
- **One-button operation:** One convenient power/program button to choose easily between 10 and 20 seconds (slow start). LED lights indicate the program change.
- **Cordless design:** The cordless design gives you maximum freedom of movement.
- **Reliable performance:** Translux Wave comes with a strong and long-lasting lithium-ion battery, with a very short recharging time. The tool-free battery exchange is quick and easy, and can be done right in the office. The built-in radiometer ensures the performance you need for the perfect restoration.

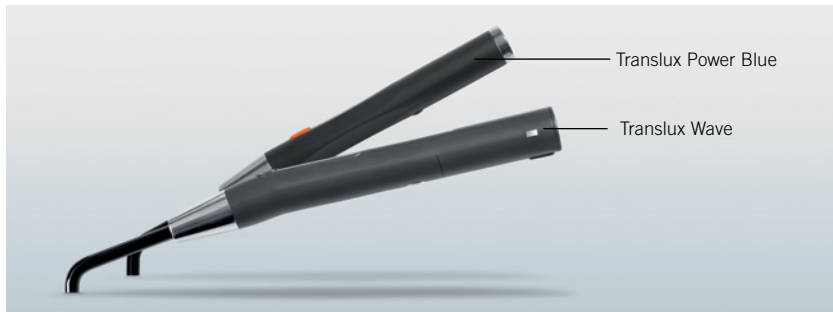

Convenience:

With its shorter and 70° angled light guide, the new Translux Wave allows easy access to molars even for patients with limited mouth-opening abilities.



Efficiency:

The battery of the handpiece and the contacts of the charging unit can easily be replaced by yourself, saving time and avoiding unnecessary down time.

Order Information

Translux Wave
Handpiece, 8 mm light guide, glare shield, charging unit, power supply
66055012

Accessories
Light guide, 8 mm, 70°
66057365

Rechargeable battery
66057363

Glare shield, orange
66016453

Replacement contacts for charger
66057382

Hand Piece incl. battery
66057383

Charging unit grey
66057384

Power supply
66057468

Technical data

Power Density
Up to 1,200 mW/cm²

Wavelength
440 – 480 nm

Dimensions charging unit
L140 x W60 x H60 mm

Dimensions hand piece
L210 mm, max. ø 23 mm

Weight
Charging unit 437 g, Hand piece: 145 g

Power Supply (charging unit)
100 to 240 V(AC), 50/60 Hz,
5 W, 12 – 18 VA

Power Supply (hand piece)
Rechargeable Lithium-ion battery

Light Source
High-luminosity mono-die LED

Light guide
ø 8 mm, 70° angle

